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Digital Camera Image Display Program
Installation／Operations Manual

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・ Installation of a camera driver

- (1) "Setup.exe" in the Driver folder of CD-ROM for installation is performed.
- (2) It chooses from "Next" in order of "install", and a driver is installed.
- (3) A USB camera will be connected if installation is completed.

Connection of a camera will start installation of a driver automatically. I carry out

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・ Installation of a digital camera image display program

- (1) "Install.exe" in CD-ROM for installation is performed.
- (2) Since an installation screen is displayed, the "Next" is chosen.
- (3) Since the screen which chooses the folder of an installation place is displayed, the "Next" is chosen.
- (4) If installation is completed, since an end screen will be displayed, "Completion" is chosen.
- (5) Installation is completion.

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
(1) The USB key is attached to a personal computer.

(2) A camera is attached to a personal computer.

(3) ASCam.exe or shortcut is chosen and performed.

(4) Starting will display the following screens.

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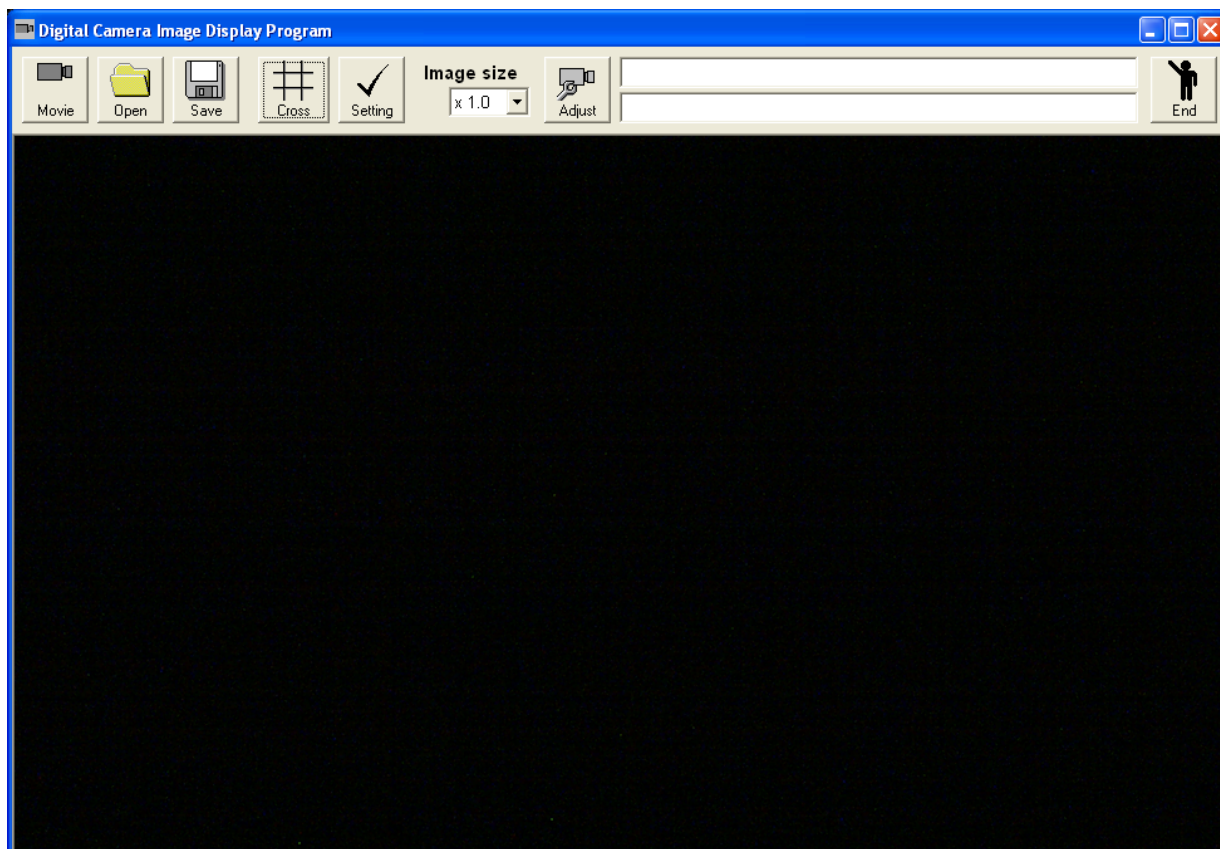※ The case where the USB key for other programs is attached and performed when it performs without attaching the USB key to a personal computer "USB key is not found.Please confirm whether a USB key is connected." and a program does not start.

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※ If the USB key is removed during execution of a program "the USB key is not found. Please

check whether the USB key is connected. It is displayed as ". If the USB key is attached and a "retry" button is clicked, a program can be used succeeding. A click of "cancellation" button will end a program. In this case, keep in mind that no image is saved.

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The "end" icon at the upper right of a screen or "x" mark is clicked.

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It is movie mode at the time of starting.

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When measuring, it is made still mode.

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
The still picture file (BMP, JPG form) saved is opened.

If a file is opened, it will change to still picture mode automatically.

The image size of a still picture file is too large to display. A check message is displayed when it cannot display.

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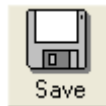A still picture (BMP, JPG form) and a movie (AVI form) are saved.

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○ The method of save of a movie

"Save" is chosen and a save place is specified.

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A file name is inputted and the kind of file is made into "a movie (*.avi)."

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If "save" is chosen, the icon of "movie" will change to "start".

"start" is pushed and save of a movie is started.

If the remaining time which chooses a "stop" for stopping or is displayed on the comment field is set to 0, it will stop automatically. (The capacity of about 2GB is consumed in 1 minute as a

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
It is an option of a setup at the time of preservation of cross hairs and a still picture.

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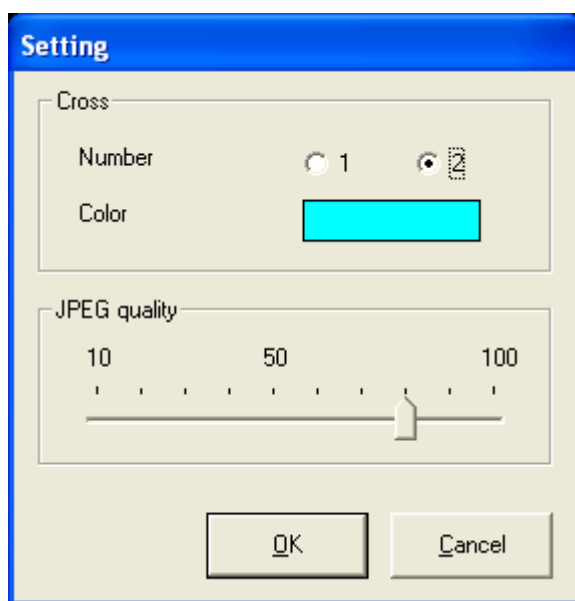
Cross hairs can perform selection of a number, and change of a color.

"JPEG quality" sets up the quality of image at the time of save.

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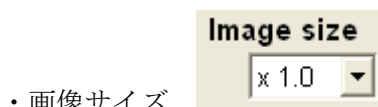It is used to indicate the image currently displayed by expansion.

The range of the magnification which can be specified is from x0.5 to x2.0.

When the image is too large, it may not indicate by expansion. Image size is set to x1.0 after a confirmation screen is displayed, when not indicating by expansion.

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A setup of "brightness", "contrast", etc., etc. can be performed at the time of movie mode.

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
A comment can be freely inputted to two lines. The inputted comment is saved as a text file of a same name at the time of still picture save. Also when image is opened, it is read simultaneously and displayed. The operation method and the remaining time are displayed at the time of save of a movie.

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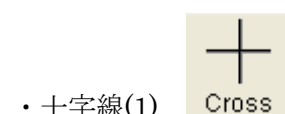
It clicks and fixes in a part to specify since cross line are displayed on a screen.

A button will be pushed in during the display of cross line.

If it clicks on an icon while on display, the display of cross line will disappear.

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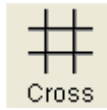Two cross line can be specified by setup.

It clicks and fixes in a part to specify since the 1st cross line are displayed on a screen. Then, since the 2nd cross line are displayed, it clicks and fixes in a part to specify similarly.

A button will be pushed in during the display of cross line.

If it clicks on an icon while on display, the display of cross line will disappear.